



# C6AHD

## Video Colposcope


---

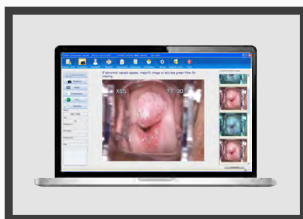

Powerful Image Display  
Data Management Software

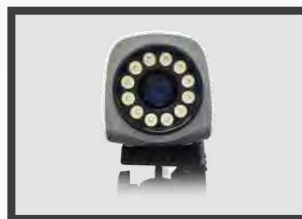

Utilizes 40 LED Lights  
for the Clearest Image

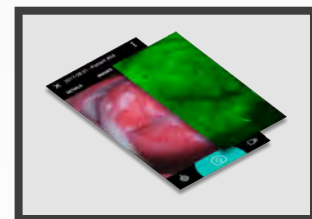

Electronic Green Filter Enhances  
the Vascular Image Display



The C6AHD Video Colposcope helps clinicians provide a comprehensive cervical examination with advanced features including a bright LED cold lighting system, high-resolution camera, and electronic green filter. With the professional software and analysis tool, patient management and reports are easy and efficient.

## Features:

- DICOM 3.0 supportive
- Real-time display of magnification (1-60x)
- Powerful Colposcopic findings gallery for comparison analysis
- LED cold lighting system reproduces an original view of the tissue
- Remote capture control and user-friendly keyboard layout for one-hand operation
- Acetic acid and iodine reaction test timer to precisely detect low/high-grade CIN
- New innovation – R-way™ Evaluation Program, which supports cervical disease diagnosis guidance for clinical users


## R-way™ Evaluation System

The R-way™ evaluation system reveals the variation of “color” on the cervix to interpret the performance of cervical diseases.

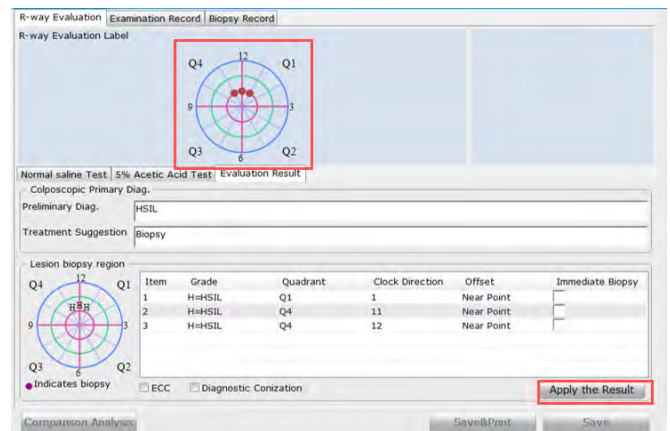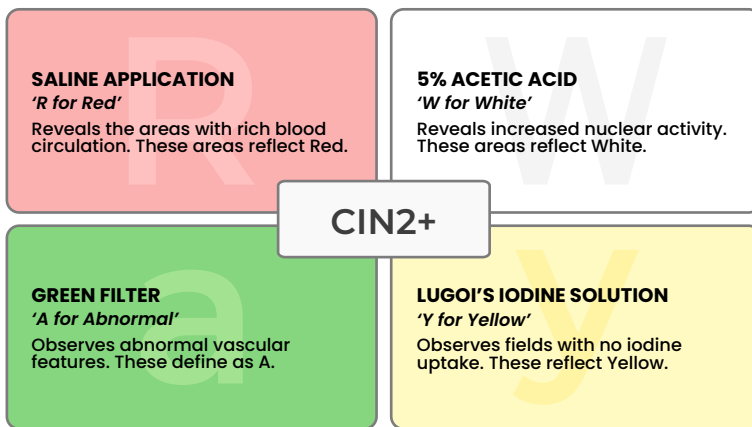

R-Way Evaluation Results

# Patient Management Software

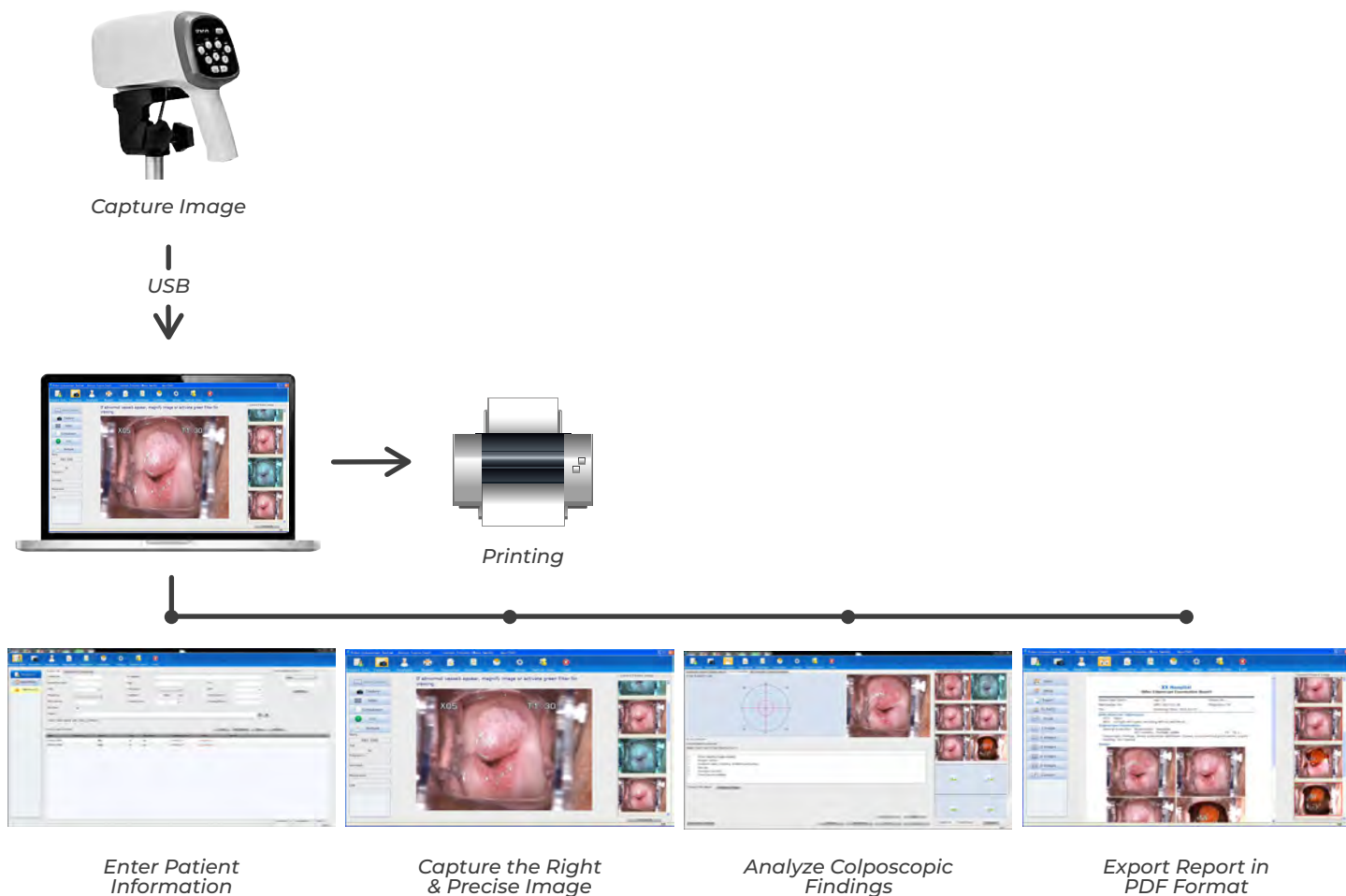

## Specifications

### DIMENSIONS

Vertical stand the chassis diameter:  $\leq \varnothing 480$  mm  
 Swing arm stand the chassis diameter:  $\leq \varnothing 550$  mm  
 MT-806 Trolley:  $\sim 500$  mm x 500 mm

### DISPLAY

Vertical stand:  $\leq 8$  kg  
 Swing arm stand:  $\leq 45$  kg (Safe Load)  
 MT-806 trolley:  $\leq 30$  kg

### TECHNICAL

Pixel: 2,130,000  
 Focus Mode: Manual, Auto  
 Output Signal: HD-Signal  
 Illumination:  $\geq 3000$  Lux  
 Magnification: 1 ~ 60X  
 Operation Distance 6X:  $\geq 200$  mm-300 mm  
 Field of View  
     3X:  $\geq \varnothing 100$  mm  
     18X:  $\geq \varnothing 15$  mm  
 System Resolution:  $\geq 900$  TVL