



The MDPro M4 is a handheld portable bladder scanner, which can connect to a phone or tablet, and can be applied in the department of surgery, emergency, ICU, obstetric, gynecology, home health, and nursing homes. It provides basis for clinical catheterization, evaluates the post-void residual after patient urination, and assists to make auxiliary diagnoses of bladder and kidney diseases.



Features

- Non-invasive 3D urinary bladder volume measurement
- Volume measurement accuracy: ±7%
- Large data storage space: 10,000+ cases
- Rapid scanning, measurement results within seconds
- No annual calibration is needed
- WiFi connection to tablet, PC, or mobile phone

Scan Modes: 3D Scan Print Mode: Optional Bluetooth Printer Working Mode(s): Easy, Expert, and Intelligence Mode

Data Transmission & Storage

Data Transmission:

- WiFi
- PC software included for PDF reports

Data Storage:

- Internal storage on phone or tablet (iOS and Android)
- 1000 patient history information review





ASC



Longterm Care




Display



Expert Mode (Displays an Ultrasonic Image)



Intelligence Mode (Displays a Real-Time Bladder Scan Projection Position)

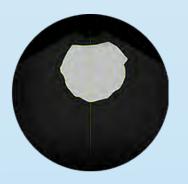

Easy Mode (Displays an Easy to View Figure)

Working Modes

Modes selection:

Users can choose between easy, expert, and intelligence mode. When performing a pre-scan, the operator will observe the following exam results:

- Move the probe until the center of the bladder image reaches the green line
- Press scan button once the biggest bladder section is detected
- Twelve real time bladder ultrasonic images will be displayed during scanning
- Volumetric reading is provided along with images and patient information for saving of test

Accessories

Highlights





Automatic Bladder Positioning



Bladder Volume 3D Scan Imaging



Over 1000 Patient History Information



Docking Station Included *Optional Carrying Bag* *Optional Cart*



Wireless Report Printing via Optional Bluetooth Printer

Specifications

Modes of Scan 3D Scan

Probe Frequency 2.5 MHZ

Volume Measurement Accuracy ±7%%, ±7ml

Scanning Time < 5 Seconds

Measure Volume 0ml-999ml

Operation Mode Easy, Expert, Intelligence Mode

Language Selection

Multiple language: English, Spanish, Swedish, Portuguese, Finnish, Chinese, French, Danish, Dutch, Norsk, German, Italian

Battery Capacity

2400 mAh

Print Mode Optional Bluetooth Printer

Data Transmission WiFi PC software included for PDF reports 1000 tests internal storage on phone or tablet (iOS and Android)

Data Storage Internal Storage on Phone or Tablet PC

Continuous Usage 2 hrs 20 min

Stand By

› 4 hrs 20 min Power Saving Mode with Automatic Shutdown Options After 5, 10, 15 or 20 Min.

